



**THUASNE**

---

## **SpryStep<sup>®</sup>** **SOLUTIONS**


Dynamic ankle-foot orthosis  
(AFO)

---


IMPROVING MOBILITY  
WITH EVERY STEP



THE TWO SHELL CONFIGURATIONS,  
ANTERIOR OR POSTERIOR, ALLOW TO ADDRESS  
DIFFERENT TYPOLOGY OF BIOMECHANICAL DEFICITS.

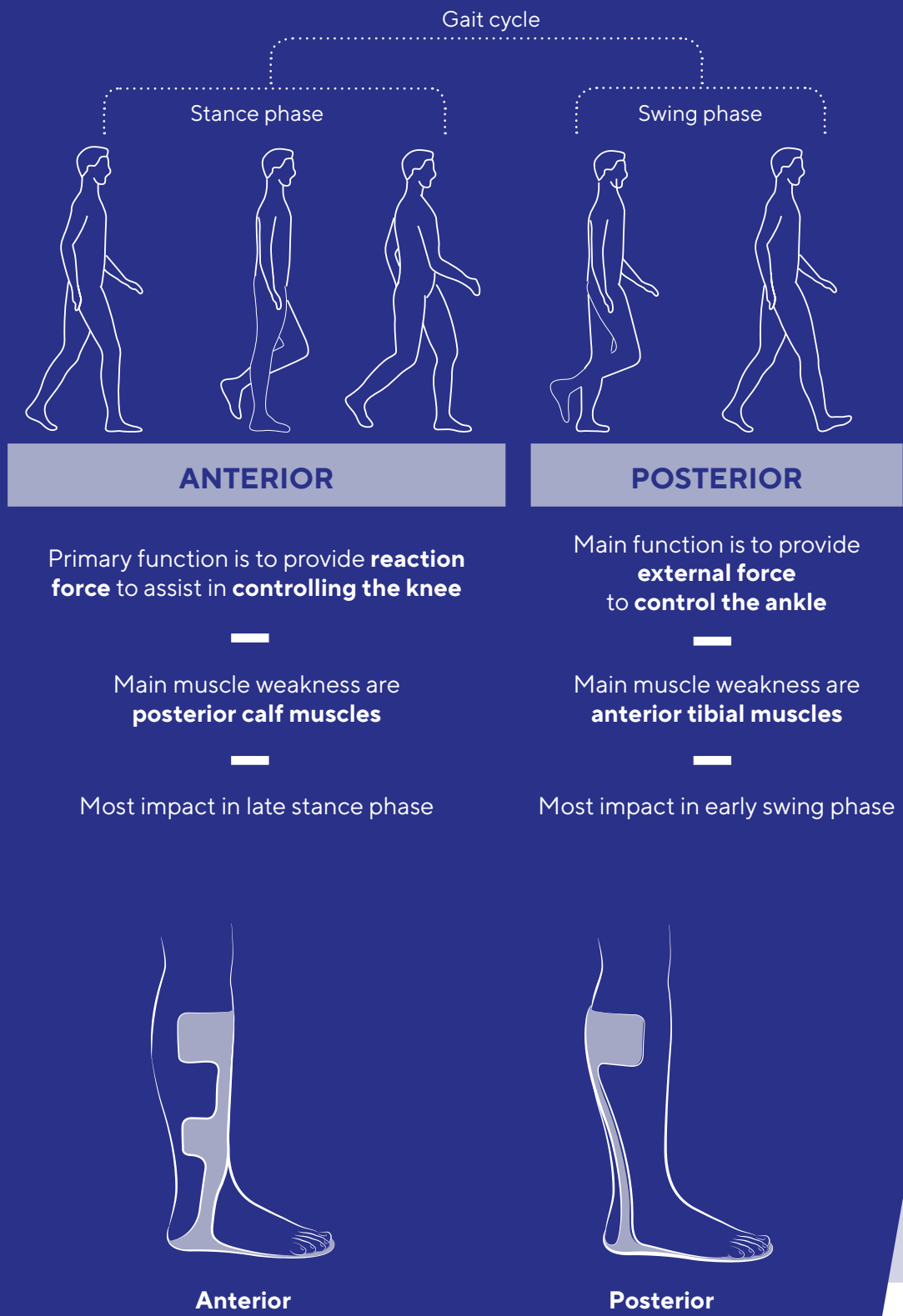

**SpryStep® A COMPLETE RANGE**  
FOR PATIENTS WITH LOWER LIMB WEAKNESS

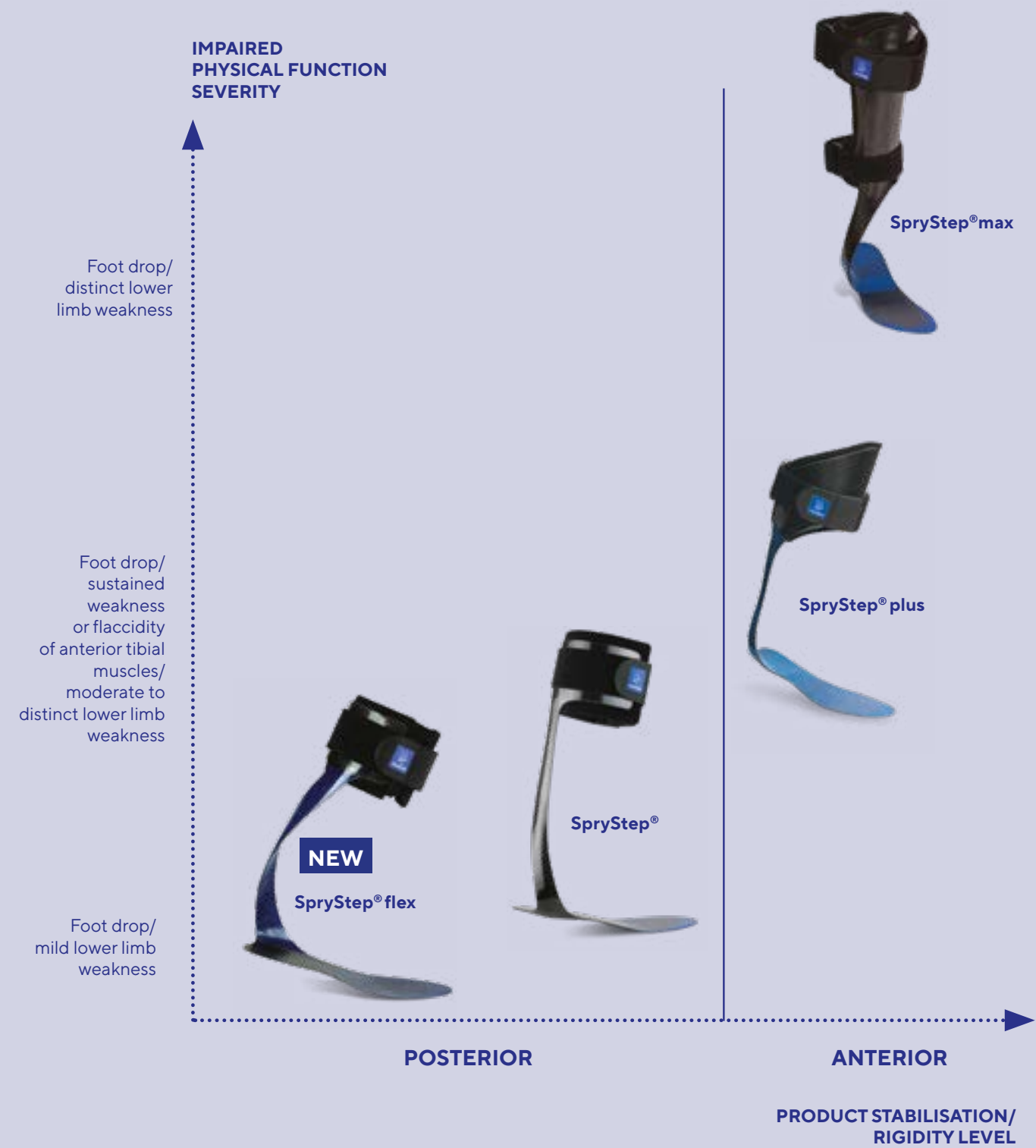

# SpryStep® RANGE:

Defined by elegance, Thuasne's SpryStep® range of AFOs is transforming the way people move. Using a proprietary blend of composite materials, the SpryStep comes in a variety of solutions that are improving mobility and quality of life for patients.

## EFFICACY

- **Energy return provided** by an optimal brace structure and the combination of high-quality composite materials.
- **Balance restoration** with a specific 3-point force system geometry which enhances foot and ankle motion and knee stabilisation
- **Improvement of gait** thanks to strut stiffness and graduated footplate flexibility

## DURABILITY

- **Durable device** that uses a well-balanced combination of composite materials in a brace structure that maximizes their properties
- **Resistant:** 2 million cycles without any compromise of the structural integrity\*  
(2 million cycles, equivalent to about 2 years of use)

## COMFORT

- **3D-shaped inner padding**, high density foam with soft surface

## PRACTITIONER FRIENDLY

- **Pre-assembled product**
- **Easy adjustment**, fitting and care: low profile device to improve compliance
  - anatomically shaped design,
  - low materials thickness
  - adapted to individual calf size thanks to removable and trimmable straps
  - washable fabric part
- **Trimmable footplate** to fit the foot length




SpryStep® flex



SpryStep®



SpryStep® plus

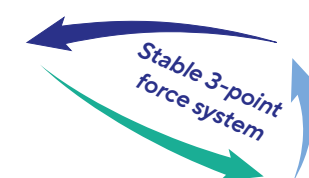


SpryStep® max



Sole pattern



Use the original sole of the shoe to determine the size of the pattern  
Trace it on the **SpryStep®** sole

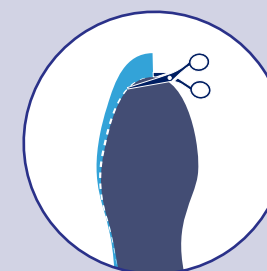

**Cutting out the SpryStep® sole**  
Use scissors to trim the **SpryStep®** sole  
**Only the blue area can be cut**

\* cycle-testing has been performed under  
ISO 10328 Servo-Pneumatic Test System





# SpryStep® flex

## INDICATIONS

### Foot drop of:

- neurological origin from central nervous system (motor impairment due to stroke, multiple sclerosis, Parkinson's disease, spinal cord injury, cerebral palsy, etc.)
- peripheral nervous system (neuropathy, nerve injuries, nerve damage, radiculopathy) or muscular origin (muscular dystrophy, congenital myopathy, etc.) without knee impairment

**Sprystep® flex is adapted for dynamic patients with high impact activities (walking, running...)**

### Posterior-lateral strut


**SpryStep® flex** is the newest member of the **SpryStep®** family. It features the same posterior-lateral strut as **SpryStep®**, but the greater spiral curve of the strut and a different material blend give the **SpryStep® flex** more flexibility.

### Design achievements and flexible properties

Adapted for the most active patients, the unique design of the **SpryStep® flex** and its flexibility provide the ability to fit patients with mild foot drop.

### Specifically designed for high-activity patients

Additionally, the most active patients (marathon runners, etc.) can benefit from its unique dynamic properties without risk of breakage or delamination.

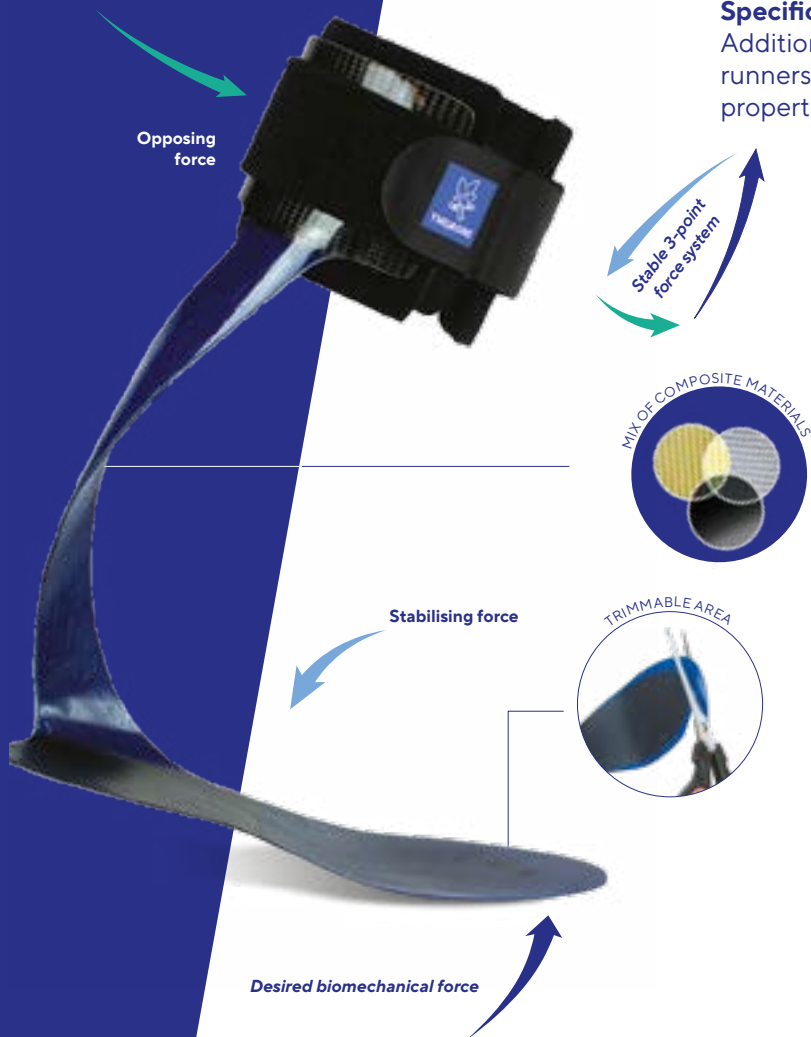

# SpryStep®

## INDICATIONS

Foot drop of neurological, traumatic or muscular origin

### Proprietary fabrication

Precise selection and layup of diverse composite materials create remarkable strength and beautiful cosmetic characteristics

### Design achievements

Geometry and material diversity achieve a necessary interaction of flexibility and stiffness during the gait cycle

### Posterior lateral strut

The spiral strut is uniquely positioned behind the malleoli to amplify energy return and fit more easily into most shoes



## SpryStep® plus



### INDICATIONS

#### Lower limb weakness of neurological, traumatic or muscular origin:

- Fatigable footdrop
- Foot slap, footdrop
- Excessive plantar flexion during swing phase (secondary to weak dorsiflexors)
- Weakness of the pretibial muscles  $\leq 3$  Plantar flexor strength 3 or 4
- Mild knee instability during stance phase
- Moderate knee instability during stance phase
- Mild quadriceps weakness
- Hip hiking

#### Design achievements


A combination of optimal geometry and well-balanced material selection delivers a “spry” gait for the patient

#### Enhanced ground reaction

The spiral strut applied to an anterior shell design provides effective dynamic reaction, especially noticeable on sloped surfaces

#### Anatomical padding

The anterior shell is lined with an elegantly engineered 3D-shaped pad featuring a recessed tibia crest relief zone

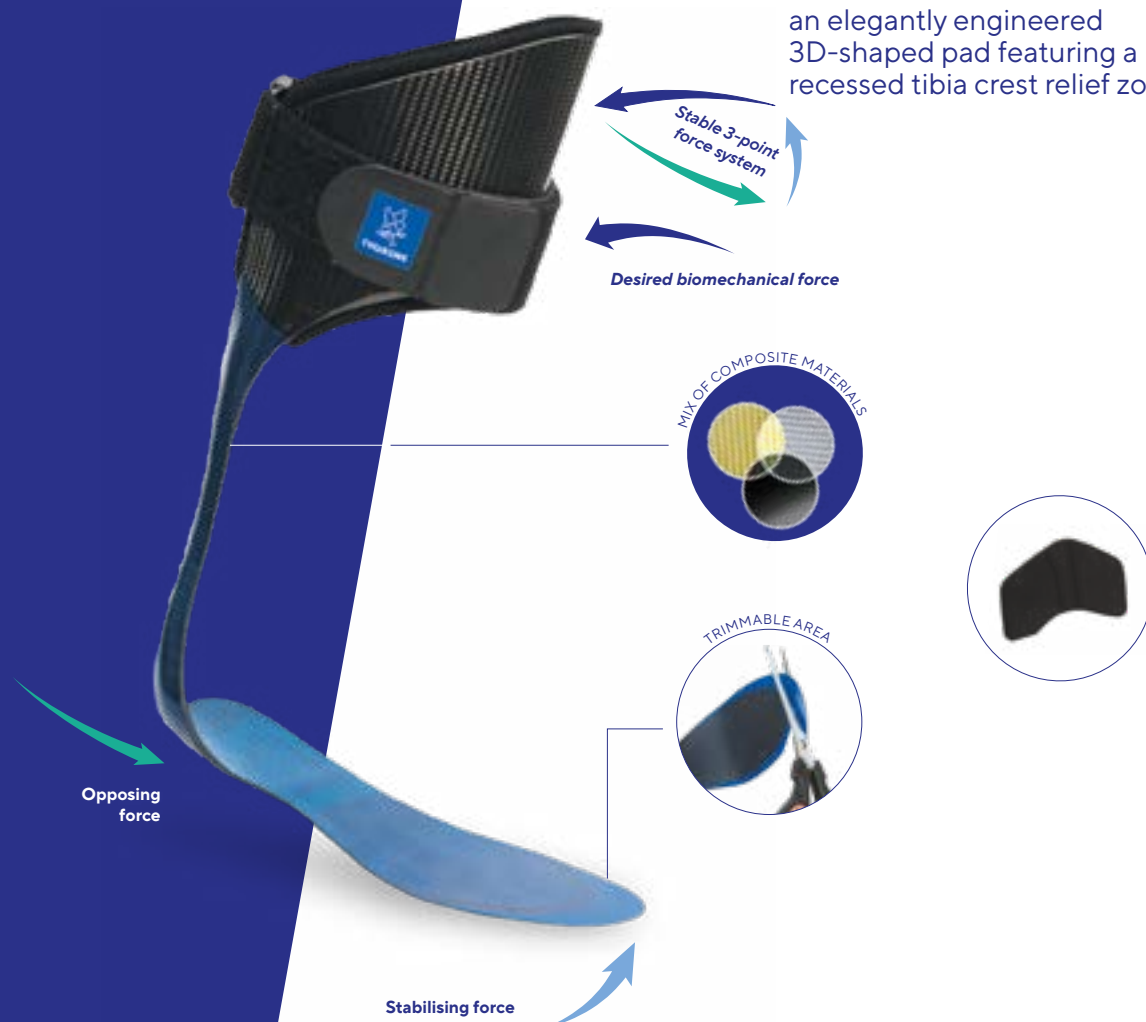

## SpryStep® max



### INDICATIONS

#### Lower limb weakness of neurological, traumatic or muscular origin:

- Excessive dorsiflexion in stance phase (due to weak plantar flexors)
- Excessive knee flexion during stance phase (secondary to weak plantar flexors)
- Quadriceps weakness
- Pain in movement towards dorsiflexion
- Foot drop
- Foot slap

#### Design achievements

The structural and material elements create graduated stiffness to achieve knee stability and enhanced forward progression

**Anterior lateral strut seamlessly integrates** into the footplate anterior to the malleoli (no sharp edge) to achieve the desired stabilisation and gait-assist

**Dynamic response and enhanced knee** stabilisation provided by the full anterior shell design







# SpryStep® AFO

## PRESCRIPTION GUIDE

The SpryStep® AFO range is designed for patients who have foot and ankle deficits. It is not suitable for patients who have active ulceration or fluctuating edema.

An elegant, durable and effective range of off-the-shelf ankle-foot orthoses. Featuring varying degrees of stiffness, posterior and anterior calf cuffs and differing strut geometries, the SpryStep® range offers an effective solution for a wide variety of patient presentations.

Resistance of AFO in Dorsiflexion

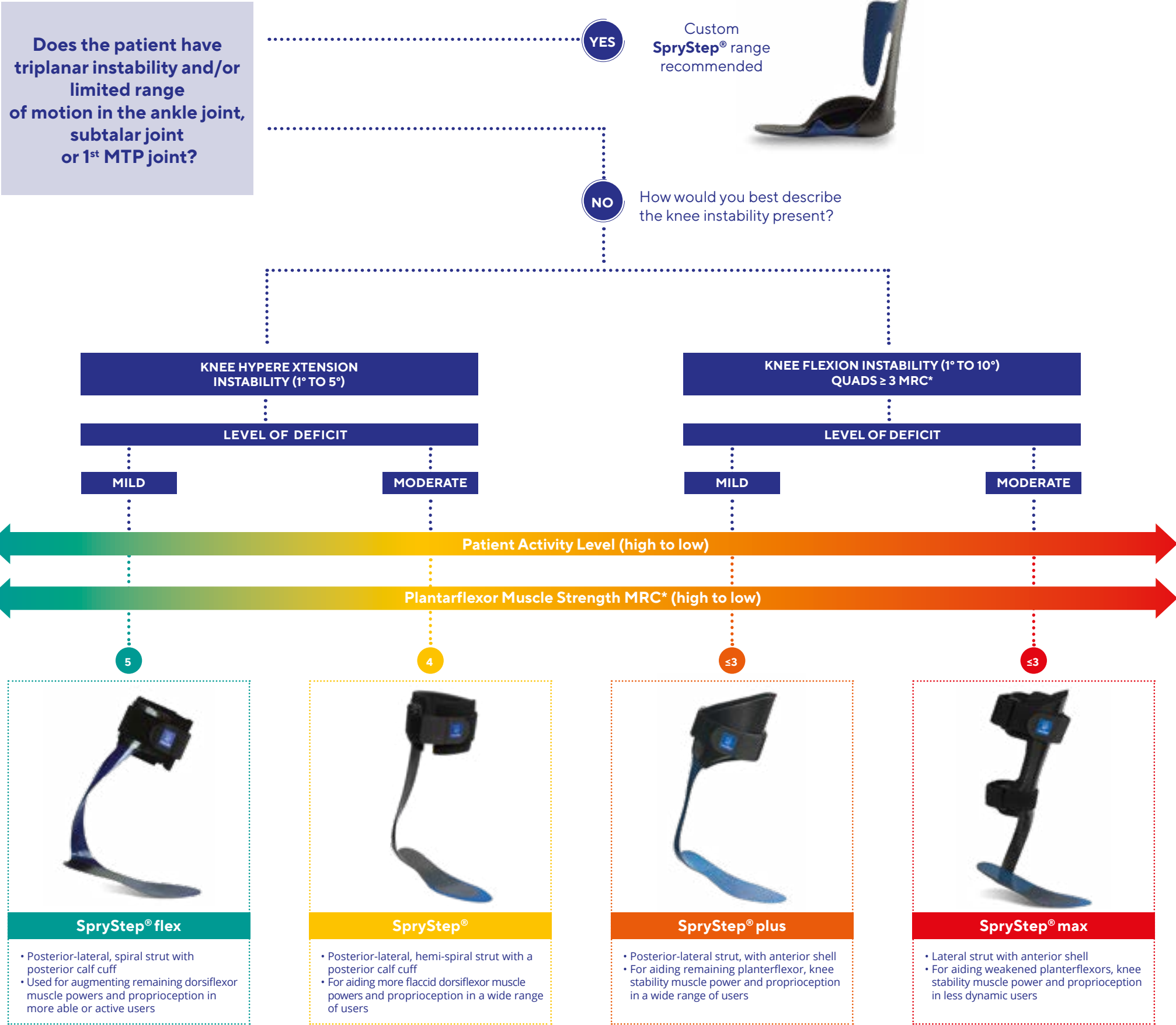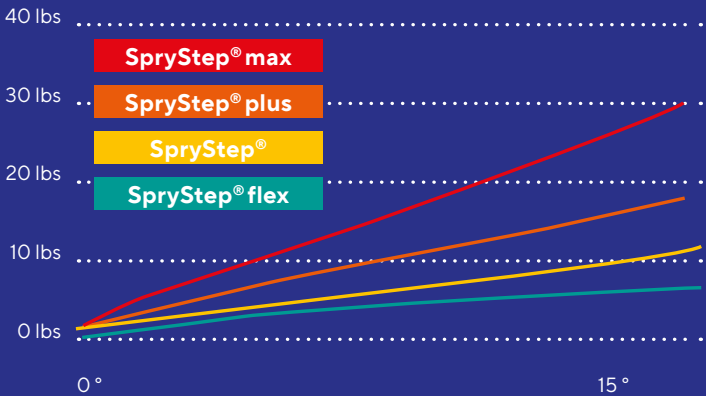

\*Medical Research Council Manual Muscle Testing Scale

# SpryStep® SOLUTIONS



## SpryStep® flex

Reference: U017 12

**2 MODELS: RIGHT AND LEFT**

| SIZE      | SHOE SIZE EU | PRODUCT HEIGHT |
|-----------|--------------|----------------|
| <b>XS</b> | 33 - 37      | 30 cm          |
| <b>S</b>  | 36 - 39      | 32 cm          |
| <b>M</b>  | 38 - 42      | 34 cm          |
| <b>L</b>  | 41 - 44      | 36 cm          |
| <b>XL</b> | 44 - 47      | 38 cm          |

For more information



## SpryStep®

Reference: U01703

**2 MODELS: RIGHT AND LEFT**

| SIZE      | SHOE SIZE EU | PRODUCT HEIGHT |
|-----------|--------------|----------------|
| <b>XS</b> | 33 - 37      | 30 cm          |
| <b>S</b>  | 36 - 39      | 32 cm          |
| <b>M</b>  | 38 - 42      | 34 cm          |
| <b>L</b>  | 41 - 44      | 36 cm          |
| <b>XL</b> | 44 - 47      | 38 cm          |

For more information



## SpryStep® plus

Reference: U017 22

**2 MODELS: RIGHT AND LEFT**

| SIZE      | SHOE SIZE EU | FOOT LENGTH |         | CALF CIRCUMFERENCE<br>25 mm BELOW FIBULA HEAD |         | PRODUCT HEIGHT |
|-----------|--------------|-------------|---------|-----------------------------------------------|---------|----------------|
|           |              | MIN         | MAX     | MIN                                           | MAX     |                |
| <b>XS</b> | 33 - 37      | 21.5 cm     | 24.5 cm | 29 cm                                         | 37 cm   | 33 cm          |
| <b>S</b>  | 36 - 39      | 23 cm       | 26 cm   | 32.5 cm                                       | 40 cm   | 36 cm          |
| <b>M</b>  | 38 - 42      | 24.5 cm     | 27.5 cm | 35.5 cm                                       | 43 cm   | 39.5 cm        |
| <b>L</b>  | 41 - 44      | 26 cm       | 29 cm   | 38.5 cm                                       | 46.5 cm | 42.5 cm        |
| <b>XL</b> | 44 - 47      | 29 cm       | 30.5 cm | 42 cm                                         | 49.5 cm | 45.5 cm        |

For more information



## SpryStep® max

Reference: U017 32

**2 MODELS: RIGHT AND LEFT**

| SIZE      | SHOE SIZE EU | FOOT LENGTH |         | PRODUCT HEIGHT |
|-----------|--------------|-------------|---------|----------------|
|           |              | MIN         | MAX     |                |
| <b>XS</b> | 31 - 34      | 19.5 cm     | 21.5 cm | 33.5 cm        |
| <b>S</b>  | 34 - 37      | 21 cm       | 23 cm   | 37 cm          |
| <b>M</b>  | 37 - 40      | 23 cm       | 25.5 cm | 39.5 cm        |
| <b>L</b>  | 40 - 43      | 25.5 cm     | 27.5 cm | 42 cm          |
| <b>XL</b> | 43 - 46      | 27.5 cm     | 29 cm   | 42 cm          |

For more information



**Help the patient to put the Spry back in its step.**

THUASNE DEUTSCHLAND GmbH  
Im Steinkamp 12 - 30938 Burgwedel - GERMANY  
www.thuasne.com

SpryStep® flex and SpryStep® are intended to provide support to the foot for patients who are suffering from foot drop of neurological, traumatic or muscular origins. SpryStep® plus and SpryStep® max are intended to provide support of the foot and the lower leg for patients who are suffering from lower limb weakness of neurological, traumatic and muscular origin. These products are regulated medical devices which bear the CE marking. Read all instructions before use and seek advice from a qualified healthcare professional. SpryStep®, SpryStep® flex, SpryStep® plus and SpryStep® max : Thuasne Deutschland GmbH, Germany  
2020 - 08 - 2003097 - © Makenzie Hollar - © Studio Caterin - © Jennifer Williams - © M-1 STUDIOS